




VASMED TECHNOLOGIES LIMITED

P.O. Box 1356, Fujairah,
Fujairah Free Zone,
Fujairah, United Arab Emirates
Tel. +971 (0) 9 222 8862
Fax. +971 (0) 9 222 8863
sales@vasmed.ae
www.vasmed.ae

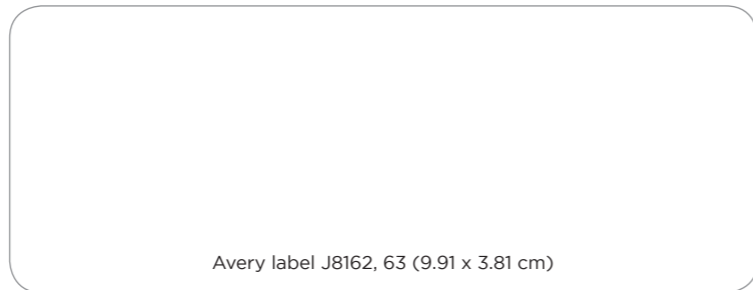

Avery label J8162, 63 (9.91 x 3.81 cm)

Torrent
Coronary and peripheral aspiration

Torrent

Aspiration Catheter



Torrent

The TORRENT Aspiration Catheter is a single operator, rapid exchange design with a dual lumen distal shaft that is compatible with 0.014” guide wire and braided proximal shaft. It has a radiopaque marker at the tip and a total working length of 140 cm. The Torrent Aspiration Catheter system comes with a preloaded removable stiffening stylet that resists kinking during delivery. The Torrent Aspiration catheter system is supplied with a luer connector connected to the aspiration line, blood filter basket, luer lock syringe and Aspiration Syringe. An appropriately sized syringe filled with an infusate may be attached to the aspiration line for the infusion of fluids.

Usable Length	140 cm
Rapid Exchange Segment	10 cms
Radiopaque marker	1 at the tip
Positioning marker	90 and 100 cm
Guide wire compatibility	0.014”
Coating	Hydro Drive hydrophilic coating

Torrent aspiration catheter designed to offer an optimal balance between crossing performance, kink resistance and thrombus-aspiration capability. Key features include :

- > Low-profile soft tip
- > Elliptical front hole aspiration
- > 25 cm long hydrophilic coated section
- > Enhanced strength from braided proximal shaft
- > Improved push ability and kink resistance
- > Removable supporting stylet.
- > Radio opaque markers at the tip for precise positioning

Compatible with 6F guiding catheter for cardiovascular use and 7F, 8F guiding catheter for peripheral use.

Ordering code – Guide catheter compatibility

Catalogue No	Version	Product Name	Diameter	Working Length	Minimum Guide Catheter Compatibility
VM-TOR-06	6F	TORRENT® Aspiration Catheter	6F	140 cm	6Fr
VM-TOR-07	7F	TORRENT® Aspiration Catheter	6F	140 cm	7Fr